



NexxZr® Shine

高效拋光。
自然光澤。

高效拋光。
自然光澤。

NexxZr® Shine



NexxZr. Shine 是一種兩步式。鑽拋光系統,用於有效拋光牙科陶瓷,例如:氧化鋯、氧化鋁、二矽酸鋰或玻璃陶瓷。

尺寸穩定、鑽石注入的橡膠拋光工具,配合簡單高效的應用得到用戶信任。此外,它們具有較長的使用壽命並為膺復體提供栩栩如生的自然光澤。NexxZr Shine 提供三種

不同結構和形狀的優化拋光工具分別是 Turn 圓盤18和, Twist 刷輪17, Tip 圓錐 5, 可使用在膺復體的整體和部分表面上。這些形狀有光滑(第 1 步,預拋光)和光澤(第 2 步,高光拋光)的顆粒度。通用的鑽石拋光膏 NexxZr Shine Paste 完善了產品系列。它可以去除微划痕並留下光滑如玻璃的表面,具有高光澤。

技術參數



類型
兩步驟鑽石
拋光系統



材料
鑽石注入橡膠拋
光工具



顆粒細度
smooth 幼細
gloss 精細



適用於
所有陶瓷



3個形狀,1個拋光膏
Turn 18, Twist 17,
Tip 5, 拋光膏

重點

1 兩步式鑽石拋光系統 適用於不同陶瓷材料,輕鬆達到油潤光澤

步驟 1
預拋光



顆粒度:
smooth 幼細



效果:
絲光



步驟 2
高光拋光



顆粒度:
gloss 精細



效果:
油潤光澤



2 耐磨混合鑽石橡膠拋光工具 栩栩如生的光澤,使用壽命長



鏡面般光滑、自然
油潤的光澤效果


3 三個協調的形狀和一個拋光膏 拋光效果佳,不受結構和表面限制



Turn圓盤 18
用於曲面



Twist 刷輪 17
用於曲面和窩溝位

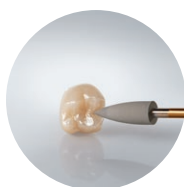

Tip圓錐 5
用於曲面和窩溝位



拋光膏
用於曲面和窩溝位



交貨形式



套裝

NexxZr Shine Starter Kit 全套裝
Turn 18 smooth / 1pcs.
Turn 18 gloss / 1pcs.
Twist 17 smooth / 2pcs.
Twist 17 gloss / 2pcs.
Tip 5 smooth / 2pcs.
Tip 5 gloss / 2pcs.
Paste 拋光膏 20 g / 1pcs.

NexxZr Shine Polishing Kit 拋光套裝
Turn 18 smooth / 1pcs.
Turn 18 gloss / 1pcs.
Twist 17 smooth / 2pcs.
Twist 17 gloss / 2pcs.
Tip 5 smooth / 2pcs.
Tip 5 gloss / 2pcs.



補充裝

NexxZr Shine Turn 18 smooth / 5pcs.
NexxZr Shine Turn 18 gloss / 5pcs.
NexxZr Shine Twist 17 smooth / 5pcs.
NexxZr Shine Twist 17 gloss / 5pcs.
NexxZr Shine Tip 5 smooth / 5pcs.
NexxZr Shine Tip 5 gloss / 5pcs.
NexxZr Shine Paste 拋光膏 20 g / 1pcs.

Manufacturer

Sagemax Bioceramics, Inc.
34210 9th Ave. South, Suite 118
Federal Way, WA 98003, USA
sagemax.com