

sagemax.com

sagemax®



NexxZr® Tool U

---

Precizno glodanje.  
Više materijala.

Precizno glodanje.  
Više materijala.

# NexxZr<sup>®</sup> Tool U



**NexxZr<sup>®</sup> Tool U su dentalni alati za glodanje visokih performansi izrađeni od vrhunskog karbida s ultra oštrim reznim rubom.**

Optimizirana geometrija reznog ruba i prostora za krhotine idealna je za obradu cirkonijeva oksida, PMMA i voska. Rezultati su iznimno glatke površine kako s tvrdim tako i s mekim CAD/CAM materijalima. NexxZr Tool U alati osvajaju korisnike svojim oštrim reznim rubovima, dobrim vijekom trajanja i atraktivnim omjerom cijene i kvalitete.

## Mehanička svojstva



**Tip**  
Radijalni rezač s  
dvostrukim rubom



**Materijal**  
Karbid vrhunske  
kvalitete



**Oštrina rubova**  
Ultra oštri



**Prikladni za**  
cirkonijev oksid,  
PMMA, vosak




**Vijek trajanja**  
Preko 120 ZrO<sub>2</sub>  
jedinica\*

## 1 Alati za glodanje od karbida vrhunske kvalitete s ultra oštrim reznim rubom


Iznimno glatke površine i dug vijek trajanja kako s tvrdim tako i s mekim CAD/CAM materijalima



**NexxZr Tool U glodalo**  
Besprijeckorni, ultraoštri  
rezní rub proizveden  
visokokvalitetnim  
postupkom brušenja

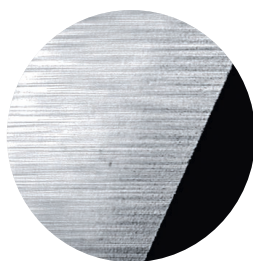

VS

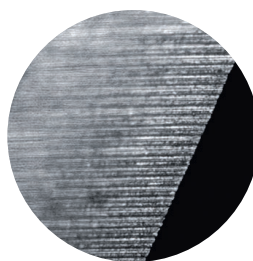

**Konvencionalno glodalo**  
Grub, neravni rezní rub  
proizveden konvencional-  
nim postupkom brušenja

## 2 Optimizirana geometrija reznog ruba i prostora za krhotine

Detaljni rezultati glodanja za  $ZrO_2$ , PMMA i vosak



CAD prikaz geometrije reznog  
ruba i prostora za krhotine



$ZrO_2$



PMMA



Wax

## 3 Atraktivan omjer cijene i kvalitete


Ekonomičan rad i visoka pouzdanost procesa

preko 120  $ZrO_2$  jedinica\*



\* Rezultati internog Sagemax ispitivanja

## Dostupno u oblicima

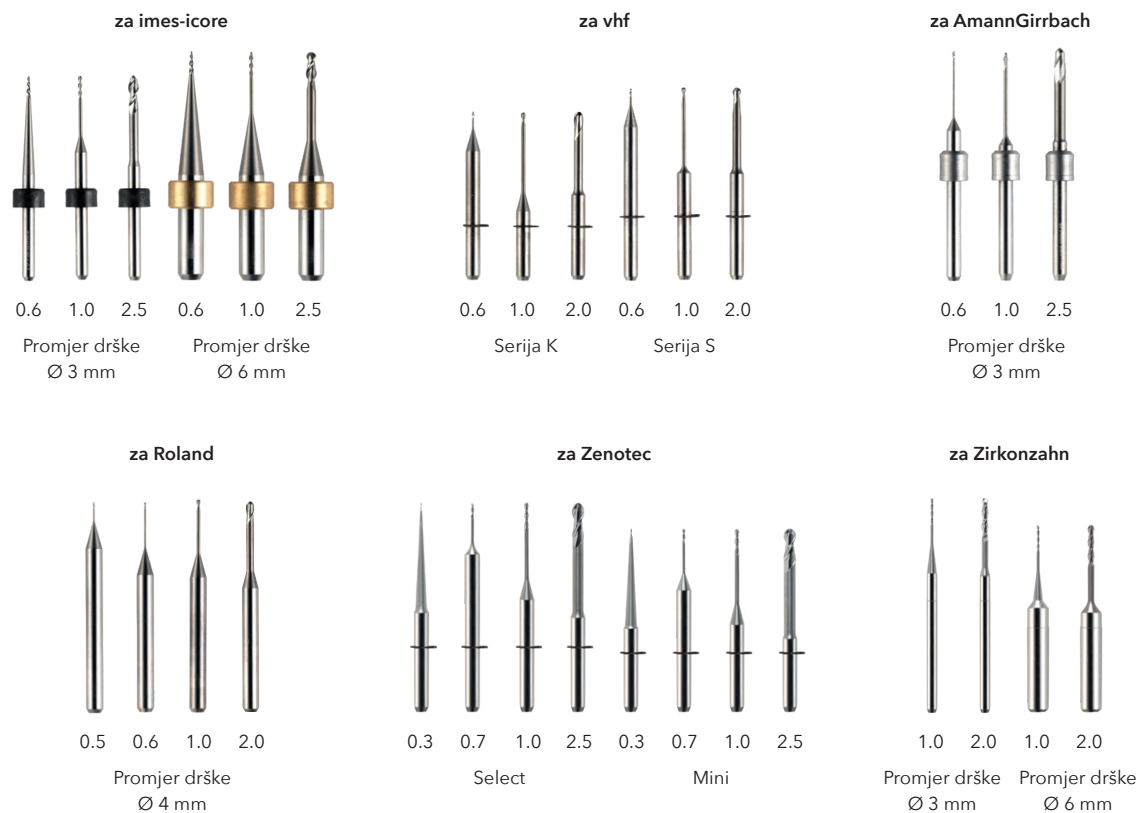


## Ostali vrhunski proizvodi



### NexxZr Tool D

Visoko kvalitetni rotirajući instrumenti za glodanje izrađeni od vrhunskog karbida s ultratvrdim, nanokristalnim dijamantnim premazom

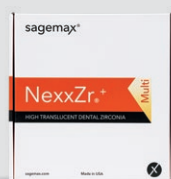

### NexxZr+ Multi

Višeslojni cirkonijev oksid s prirodnom estetikom



### PMMA Disc Multi

Visoko umreženi polimetil metakrilat s integriranim gradijentom boje za CAD/CAM-om izrađene privremene radove



### Wax Disc

Sintetički voštani disk za CAD/CAM izradu voštanih predmeta

### Manufacturer

**Sagemax Bioceramics, Inc.**

34210 9th Ave. South, Suite 118

Federal Way, WA 98003, USA

E: [info@sagemax.com](mailto:info@sagemax.com)

[sagemax.com](http://sagemax.com)